

Infyna

Ready to Use Hydrophilic Intermittent Catheter

Using a catheter
shouldn't make
people feel
like patients.



Continence Care
People First

Infyna

Ready to Use Hydrophilic Intermittent Catheter

Change the
experience
with Infyna.

Ready. Easy.




Continence Care
People First

Redefine what it feels like for your patients to use a catheter with Infyna

We listened to you and created a catheter that helps meet the needs of your patients.




Here are a few ways we redesigned the experience to help your patients feel less like patients:

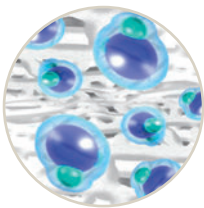

Ready to use with reduced spill¹

- The Vaporphilic technology provides a fully pre-hydrated, ready to use hydrophilic catheter, with no need to add water
- The vapor strip releases sterile water vapor that activates the catheter coating for lubrication. Once opened, the product is pre-hydrated and ready to use.

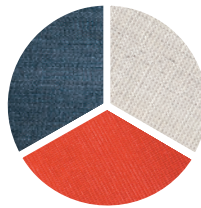

Leaves no stain

- On common clothing fabrics, specifically denim, linen, and polyester²



No-touch catheterization

- Designed for more control and smooth insertion



Open with confidence

- Easy-to-open package
- Easy tear strip opening

Available in a range
of options for women

including standard and closed systems

Infyna Ready to Use Hydrophilic Intermittent Catheter

Size (Charriere)	17 cm Nelaton	17 cm Nelaton Closed System	Quantity per Box
8 Ch	88081	86081	30
10 Ch	88101	86101	30
12 Ch	88121	86121	30
14 Ch	88141	86141	30



This intermittent catheter is a flexible tubular device that is inserted through the urethra by female and pediatric patients who need to drain urine from the bladder.

- Ready to use, single-use intermittent catheter
- Hydrophilic-coated catheter not made with phthalates
- 2 smooth catheter eyelets
- Color-coded funnel
- Not made with natural rubber latex

WARNING: To help reduce the potential for infection and/or other complications, do not reuse.

WARNING: If discomfort or any sign of trauma occurs, discontinue use immediately and consult your healthcare practitioner.

PRECAUTION: Please consult your healthcare practitioner before using this product if any of the following conditions are present:

- Severed urethra • Unexplained urethral bleeding • Pronounced stricture • False passage • Urethritis – inflammation of the urethra

PRECAUTION: Self-catheterization should follow the plan of care and advice given by your healthcare practitioner and be carried out only in accordance with the instructions for use provided. Because catheterization frequency varies by person, the recommended frequency of your catheterization should be provided by your healthcare practitioner. For any other questions about your catheterization, please contact your healthcare practitioner.

NOTE: Store boxes in a flat position and at normal room temperature.

This product consists of: polyvinylchloride funnel and tubing not made with phthalates.

REFERENCES

1. Data on file. Loose water in Onli and competitor products packaging study. Hollister. 2015.
2. Data on file. Analysis of potential loose liquid to stain clothing for Onli and competitor products. Hollister. 2015.